

BLD-015

REPORTED June 28, 2026

CLIENT

NAME Peptide Test Canada
CONTACT PTC Operations
EMAIL info@peptidetestcanada.ca
ADDRESS Canada

SCAN TO VERIFY



PTC-NEFQYXGG

END CUSTOMER

NAME Jinlong Maple International
EMAIL —
PHONE —
ADDRESS —

§ 01 Sample Information

FOR RESEARCH PURPOSES ONLY

Sample ID:	BLD-015	Total Mass:	124.8 mg	CAS:	—
Compound:	KLOW	Appearance:	Lyophilized powder	Formula:	—
Report No.:	COA-2026-SC-00120	Received:	06/24/2026	Mol Weight:	—
Lot Number:	COA-2026-SC-00120	Analyzed:	06/28/2026		

§ 02 Analytical Results

Qualitative and Quantitative chemical analysis by Liquid Chromatography Tandem Mass Spectrometry (LC-MS/MS)


PEPTIDE	EXPECTED	MEASURED		PURITY	STATUS
KPV	10.00 mg	10.53 mg	—	98.9 %	CONFIRMED
GHK-Cu	50.00 mg	46.77 mg	—	98.5 %	CONFIRMED
BPC-157	10.00 mg	9.90 mg	—	99.0 %	CONFIRMED
TB-500	10.00 mg	9.62 mg	—	98.6 %	CONFIRMED

PRECURSOR [M+H]⁺

343.20 Da

Quantifier transition ● MRM

343.20 > 226.20 CE 15 V

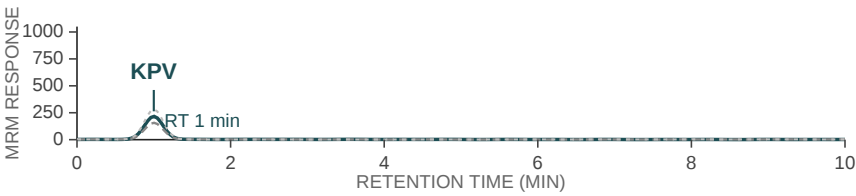

ROLE	TRANSITION	CE
● Quantifier	343.2 > 226.2	15 V
■ Qualifier	343.2 > 70.1	22 V


Agilent Poroshell 120 EC-C18 · 2.1×100 mm · 2.7 μm · Solvent A: 0.1 % formic acid in H₂O · Solvent B: 0.1 % formic acid in Acetonitrile

PRECURSOR [M+H]⁺

341.20 Da

Quantifier transition ● MRM

341.20 > 156.10 CE 15 V

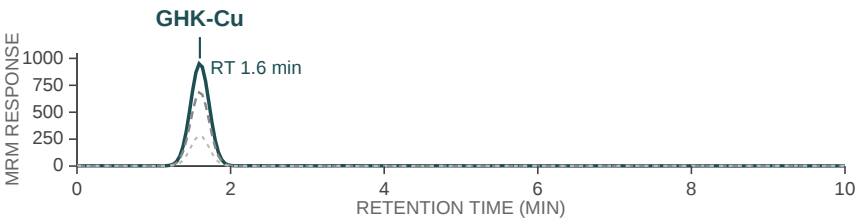

ROLE	TRANSITION	CE
● Quantifier	341.2 > 156.1	15 V
■ Qualifier	341.2 > 110.1	25 V


Agilent Poroshell 120 EC-C18 · 2.1×100 mm · 2.7 μm · Solvent A: 0.1 % formic acid in H₂O · Solvent B: 0.1 % formic acid in Acetonitrile

PRECURSOR [M+2H]²⁺

710.40 Da

Quantifier transition ● MRM

710.40 > 70.10 CE 22 V

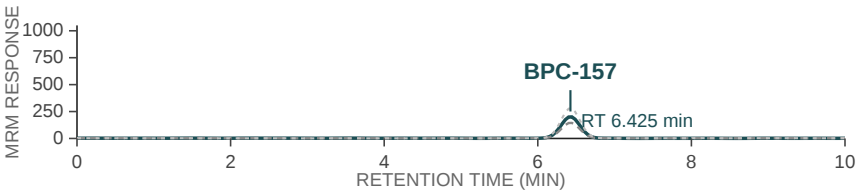

ROLE	TRANSITION	CE
● Quantifier	710.4 > 70.1	22 V
■ Qualifier	710.4 > 251.1	18 V


Agilent Poroshell 120 EC-C18 · 2.1×100 mm · 2.7 μm · Solvent A: 0.1 % formic acid in H₂O · Solvent B: 0.1 % formic acid in Acetonitrile

PRECURSOR [M+2H]²⁺

445.30 Da

Quantifier transition ● MRM

445.30 > 369.30 CE 12 V

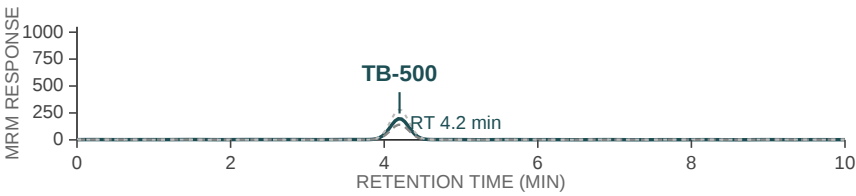

ROLE	TRANSITION	CE
● Quantifier	445.3 > 369.3	12 V
■ Qualifier	445.3 > 241.2	15 V

Agilent Poroshell 120 EC-C18 · 2.1×100 mm · 2.7 μm · Solvent A: 0.1 % formic acid in H₂O · Solvent B: 0.1 % formic acid in Acetonitrile

DISCLAIMER · RESEARCH USE ONLY

This report documents **analytical results only**. Provided for **research purposes only**; not a drug-release or quality certification. SideChain Analytics makes no representation regarding resale, clinical use, or human / veterinary administration of the material.